

KRIBIOLISA™ *E.coli* Host Cell DNA Kit






REF : KBBR01

Pack Size: 96 wells

Ver 5.0


RUO

DNA Assay Kit for Accurate Quantitation of *E.coli* Host Cell DNA from cell culture supernatant and other biological preparations

RUO	For Research Use Only	REF	Catalog Number
	Store At	LOT	Batch Code
	Manufactured By		Biological Risk
	Expiry Date		Consult Operating Instructions

For Research Purposes Only. Purchase does not include or carry the right to resell or transfer this product either as a stand-alone product or as a component of another product. Any use of this product other than the permitted use without the express written authorization of KRISHGEN BioSystems is strictly prohibited.

REF KBBR01

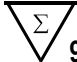 96 wells



KRISHGEN BioSystems | For US / Europe: toll free +1(888)-970-0827 tel: +1(562)-568-5005
For Asia / India: tel: +91(22)-49198700
Email: sales1@krishgen.biz

KRIBIOLISA™ *E.coli* Host Cell DNA Kit

Introduction:

For the production of biotechnological products, expression of therapeutic proteins in *E.coli* cells is commonly used. Due to which there is a possibility of DNA contamination from the host cells in the products which are manufactured. However, the limit of DNA contamination from the cell lines has been set by regulatory agencies. According to WHO the limit of DNA contamination should not exceed 10 ng/dose in the products manufactured. One of the most common method which is widely used for detection of host cell DNA contamination by biopharmaceutical manufacturers is the PicoGreen™ Dye based assay. The research and development department at Krishgen Biosystems has designed a kit to detect and quantify host cell DNA in the products manufactured by recombinant expression in *E.coli* cell lines based on PicoGreen™ Dye.

Principle:

The Krishgen *E.coli* Host Cell DNA Kit is based on the DNA dye binding assay which utilizes PicoGreen™ dye. The DNA samples and the standards are reacted with PicoGreen™ dye, which is a DNA intercalator that binds strongly to the double stranded DNA. Upon binding to DNA, PicoGreen™ dye fluoresces with an excitation of 485 nm and emission of 525 nm. The intensity of the fluorescent signal is proportional to the quantity of DNA in the standard or the samples.

Materials Provided:

DNA Extraction Reagents	
*Buffer AL (Lysis solution)	35 ml
Pico Green Assay	
<i>E. Coli</i> DNA Standards 20 ug/ml	100 ul
PicoGreen™ dsDNA reagent	210 ul
TE Buffer (10X)	5 ml
Assay Plate with Plate Sealing Foil	1 x 96 wells
Instruction Manual	1

****Buffer AL may precipitate at cold temperature. If so, dissolve it in 37°C water bath.***

Materials Required By the End User:

Reagent

1. Isoamyl Alcohol
2. Chloroform
3. Phenol
4. 70% ethanol.
5. Isopropanol
6. Optional RNase solution (not provided)

Disposable Material

1. Sterile 1.5 ml microcentrifuge tubes
2. RNase free pipette tips

Equipment

1. Micro centrifuge (4°C)
2. Adjustable and properly calibrated pipettes to measure volumes ranging from 10ul to 1000ul.
3. Vortex
4. Heating Block; 37°C and 65°C
5. Appropriate personal protector (e.g. lab coat, disposable gloves, goggles etc.)
6. Fluorimeter

KRIBIOLISA™ *E.coli* Host Cell DNA Kit

Handling/Storage

1. Buffer AL should be stored at room temperature (15-25°C). A precipitate can be formed in Buffer AL under cool condition. In such case, heat the bottle at 56°C until the precipitate completely dissolves.
2. Buffer PP should be stored at 2-8°C and bring to Room Temperature when in use.
3. The PicoGreen™ dsDNA reagent is light sensitive and should be stored in the dark until use at -20 °C.
4. *E.coli* DNA standards should be stored at -20 °C.
5. All the components in the kit are stable until the expiration date mentioned in the kit.

Health Hazard Warnings

1. All the reagents provided may be harmful if ingested, inhaled or absorbed through the skin.
2. For Research or Manufacturing Use only.

Precautions

1. The DNA extraction reagents, Buffer AL and PP contain irritant which is harmful when it comes in contact with skin or eyes or when inhaled or swallowed Care should be taken during handling. Always wear gloves and eye protector, and follow standard safety precautions.
2. The PicoGreen™ Solution is a DNA intercalating dye and should be handled with care. When working with PicoGreen™ Solution, use gloves, eye protector and mask to avoid contact with skin or clothing and inhalation of vapor. In case of contact, wash immediately with plenty of water and seek medical advice.

Procedural Notes

1. Bring the PicoGreen™ solution to room temperature before proceeding with the experiment.
2. The PicoGreen™ Solution will bind to all double stranded DNA. It is important to keep the working area clean to avoid contamination by DNA in the environment. Thoroughly clean pipettes and the immediate working area before initiating the protocol.
3. Always make sure the centrifuge is balanced to ensure proper assay performance. An improperly balanced centrifuge can result in loose pellets, which can adversely affect recovery and assay precision

Reagent Preparation (all reagents should be diluted immediately prior to use):

1. **TE buffer (1X)**
To make TE (1X) buffer, add 2 ml of TE (10X) buffer to 18 ml of double distilled water.
2. **Preparation of PicoGreen™ Solution**
Dilute the PicoGreen™ solution in the ratio of 1:100. For example, if the experiment has 20 samples, add 50 ul of PicoGreen™ to 4950 ul of pre-warmed (37°C) TE Buffer. Mix by gently vortex for 5 seconds.
Note: PicoGreen™ Solution must be diluted fresh for each experiment.
3. **Preparation of PP Buffer**
PP buffer is prepared by adding highly pure Phenol: Chloroform: Isoamyl Alcohol in the ratio of 25:24:1 and the pH is adjusted to 8.05- 8.35 and should be stored at amber bottles.
Note: At 15 -30°C the buffer will appear as a one clear solution (homogeneous phase).

Sample Preparation: Protocol for Extraction of DNA

1. Dilute all test samples to concentrations within the analytical range of the assay and to < 20 mg/ml total protein using TE buffer (1X). All samples should be diluted at least 1:2.
2. Transfer 300 ul of Buffer AL to a fresh 1.5 ml micro centrifuge tube.
3. Add 300 ul of sample to the tube containing buffer AL and pipette 5 to 6 times to resuspend thoroughly.
4. Incubate the lysate at 37°C until clumps disappear.

KRIBIOLISA™ *E.coli* Host Cell DNA Kit

Note: Generally, cell lysis is completed in 5 min. However, complete resuspending of the cells is crucial for good DNA yield. If the clumps are still visible after 1 hour, add an additional 100ul of Buffer AL and repeat incubation.

5. If RNA-free DNA is required, add 1 ul of RNase solution (20 mg/ml) to the lysate and mix the sample by inverting the tube 5 times. Incubate the mixture for 15 min at 37°C.
6. Cool the sample to room temperature. Add 100 ul of Buffer PP to the mixture and vortex vigorously for 15sec. Centrifuge for 2 min at 14,000 x g.
Note: Incubate the sample on ice for 5min before centrifugation. This may increase the yield of DNA. A dark brown protein pellet should be visible.
7. Carefully transfer the supernatant to a fresh 1.5 ml micro centrifuge tube containing 300 ul of isopropanol (room temperature) and gently mix the solution by inversion until the white thread-like strands of DNA forms a visible mass.
Note: Do not vortex after addition of isopropanol.
8. Centrifuge at 14,000 x g for 1 min. Decant the supernatant and add 300 ul of 70% ethanol (room temperature). Gently invert the tube several times to wash the DNA pellet and the side walls of the tube.
9. Centrifuge at 14,000 x g for 1 min. carefully discard the ethanol by aspirating or pipetting. Invert the tube on clean absorbent paper and air dry the pellet for 10 to 15 min.
Note: The DNA pellet is very loose at this point and care must be taken to avoid missing the pellet. Ethanol should be completely removed, but over drying will make the rehydration of DNA pellet difficult.

Assay Procedure:

Detection and Quantification of *E.coli* Host Cell DNA

1. Bring all reagents to room temperature prior to use. It is strongly recommended that all standards and samples be run in duplicate or triplicate. A standard curve is required for each assay.
2. **Standards Preparation:** Dilute 5 ul of original Standard (20 ug/ml) with 495 ul of TE buffer generate a 200 ng/ml Standard stock solution. Thus, the *E.coli* DNA standard concentrations are 200 ng/ml, 100 ng/ml, 50 ng/ml, 25 ng/ml, 12.5 ng/ml, 6.25 ng/ml, 3.125 ng/ml, 1.56 ng/ml and TE buffer (1X) serves as the zero standard (0ng/ml). Prepare the Standards by serially diluting the standard stock solution as per the below table.

Standard Concentration	Standard Vial	Dilution Particulars
200 ng/ml	Standard No.8	5 ul Standard Provided (20 ug/ml) + 495 ul TE buffer (1X)
100 ng/ml	Standard No.7	250 ul Standard No.8 + 250 ul TE buffer (1X)
50 ng/ml	Standard No.6	250 ul Standard No.7 + 250 ul TE buffer (1X)
25 ng/ml	Standard No.5	250 ul Standard No.6 + 250 ul TE buffer (1X)
12.5 ng/ml	Standard No.4	250 ul Standard No.5 + 250 ul TE buffer (1X)
6.25 ng/ml	Standard No.3	250 ul Standard No.4 + 250 ul TE buffer (1X)
3.125 ng/ml	Standard No.2	250 ul Standard No.3 + 250 ul TE buffer (1X)
1.56 ng/ml	Standard No.1	250 ul Standard No.2 + 250 ul TE buffer (1X)

3. Add 250 ul of Standards and Samples to the micro centrifuge tube.
Note: If you want to check sensitivity on your end you can dilute samples further. Please check our justification summary
4. Add PicoGreen™ in 1:1 proportion to standards. E.g. Add 250 ul of pre-diluted and pre-warmed PicoGreen™ reagent to the 250 ul of Standards. For samples directly add 250 ul of PicoGreen™ reagent in the DNA Pellet obtained and mix to dissolve the pellet. Remember to run your each sample and standards in duplicate.
5. Vortex each tube for 5 seconds and incubate for 5 minutes at Room Temperature under dark.
6. Remove 100 ul from each tube and pipette into the Assay Plate (provided). Seal the plate with the foil provided and incubate for 5 minutes at room temperature.

KRIBIOLISA™ *E.coli* Host Cell DNA Kit

7. Read the plate at Excitation wavelength: 485nm, Emission wavelength: 525nm, Cutoff: 515nm

Calculation of Results:

The standards may be used to construct a standard curve with values reported in ng/ml host cell DNA. This data reduction may be performed through computer methods using curve-fitting routines such as point-to-

point, Cubic spline and 4 parameter logistic fit. Do not use linear regression analysis to interpolate values for samples as this may lead to significant inaccuracies.


Note:

1. Precision on triplicate samples should yield average % coefficients of variation of less than 15% for samples in the range of 500pg/ml to 150ng/ml.
2. Coefficients of variation for samples above or below this range (500pg/ml to 150ng/ml) may be greater than 15%.

Typical Data:

E.coli DNA Standard Concentration (ng/ml)	Relative Fluorescence values (RFU)	Relative Fluorescence values after Blanking (RFU)
0	104	0
1.56	227	123
3.156	344	240
6.25	683	579
12.5	981	877
25	2253	2149
50	4584	4480
100	6336	6232
200	15723	15619

Typical Graph

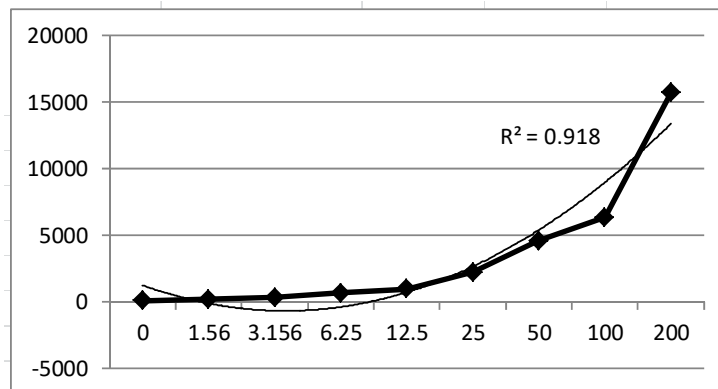

Quality Control:

E.coli Host Cell DNA Kit is manufactured in strict clean condition and its degree of cleanness for the components provided is monitored periodically. For quality control, the quality certification process is carried out thoroughly and only the qualified is delivered.

It is recommended that for each laboratory assay appropriate quality control samples in each run to be used to ensure that all reagents and procedures are correct.

Precautions:


Do not mix reagents from different kits or lots. Reagents and/or antibodies from different manufacturers should not be used with this set.

KRIBIOLISA™ *E.coli* Host Cell DNA Kit


Specificity of the Picogreen:

Fluorescent enhancement of the Picogreen Reagent on binding to dsDNA Standards.

The Ecoli DNA standards used in the kit have been validated for binding and sensitivity to Picogreen dye. The assay is linear and validated over three orders of magnitude permitting an accurate measurement of the DNA from Ecoli Host Cells.

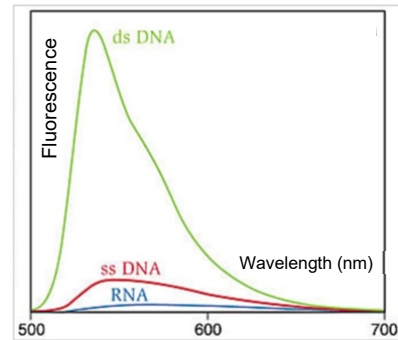

Fluorescence spectra:

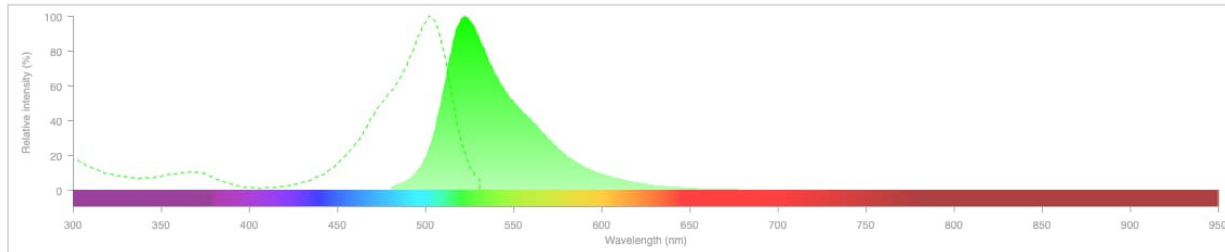

Precision:

DNA samples were prepared in a Human IgG matrix at various concentrations spanning the Standard Curve. Four preparations were made for each sample and duplicate wells were collected for each preparation.

Pool	Concentration	% Nominal	Intra-Assay CV%	Inter-Assay CV%
Low	1 ng/ml	98%	20.1%	14.0%
Medium	10 ng/ml	95%	12.6%	10%
High	200 ng/ml	92%	10%	9.7%

Safety Precautions:

- This kit is **For Research Purposes only**. Follow the working instructions carefully.
- The expiration dates stated on the kit are to be observed. The same relates to the stability stated for reagents
- Do not use or mix reagents from different lots.
- Do not use reagents from other manufacturers.
- Avoid time shift during pipetting of reagents.
- All reagents should be kept in the original shipping container.
- Some of the reagents contain small amount of sodium azide (< 0.1 % w/w) as preservative. They must not be swallowed or allowed to come into contact with skin or mucosa.
- Source materials maybe derived from human body fluids or organs used in the preparation of this kit were tested and found negative for HBsAg and HIV as well as for HCV antibodies. However, no known test
- Guarantees the absence of such viral agents. Therefore, handle all components and all patient samples as if potentially hazardous.
- Since the kit contains potentially hazardous materials, the following precautions should be observed
 - Do not smoke, eat or drink while handling kit material
 - Always use protective gloves
 - Never pipette material by mouth
 - Wipe up spills promptly, washing the affected surface thoroughly with a decontaminant.
- In any case GLP should be applied with all general and individual regulations to the use of this kit.



KRIBIOLISA™ *E.coli* Host Cell DNA Kit

LIMITED WARRANTY

Krishgen Biosystems does not warrant against damages or defects arising in shipping or handling, or out of accident or improper or abnormal use of the Products; against defects in products or components not manufactured by Krishgen Biosystems, or against damages resulting from such non-Krishgen Biosystems made products or components. Krishgen Biosystems passes on to customer the warranty it received (if any) from the maker thereof of such non Krishgen made products or components. This warranty also does not apply to Products to which changes or modifications have been made or attempted by persons other than pursuant to written authorization by Krishgen Biosystems.

THIS WARRANTY IS EXCLUSIVE. The sole and exclusive obligation of Krishgen Biosystems shall be to repair or replace the defective Products in the manner and for the period provided above. Krishgen Biosystems shall not have any other obligation with respect to the Products or any part thereof, whether based on contract, tort, and strict liability or otherwise. Under no circumstances, whether based on this Limited Warranty or otherwise, shall Krishgen Biosystems be liable for incidental, special, or consequential damages.

This Limited Warranty states the entire obligation of Krishgen Biosystems with respect to the Products. If any part of this Limited Warranty is determined to be void or illegal, the remainder shall remain in full force and effect.

Krishgen Biosystems. 2024

THANK YOU FOR USING KRISHGEN PRODUCT!

KRISHGEN BIOSYSTEMS®, GENLISA®, DHARMAPLEX™, GENBULK™, GENLISA™, KRISHZYME®, KRISHGEN®, KRIBIOLISA®, KRISHPLEX®, TITANIUM®, QUALICHEK® are registered trademarks of KRISHGEN BIOSYSTEMS. ©KRISHGEN BIOSYSTEMS. ALL RIGHTS RESERVED.

KRISHGEN BIOSYSTEMS | OUR REAGENTS | YOUR RESEARCH |