
Human NT-Terminal Fragment of The Medullary Natriuretic Peptide, NT-ProBNP GENLISA™ ELISA

: KBH4693 REF Ver 3.0 **RUO**

Enzyme Immunoassay for the Quantitative Determination of NT-Terminal Fragment of The Medullary Natriuretic Peptide, NT-ProBNP in Human serum, plasma and other biological samples.



For Research Use Only. Purchase does not include or carry the right to resell or transfer this product either as a stand-alone product or as a component of another product. Any use of this product other than the permitted use without the express written authorization of KRISHGEN BioSystems is strictly prohibited.

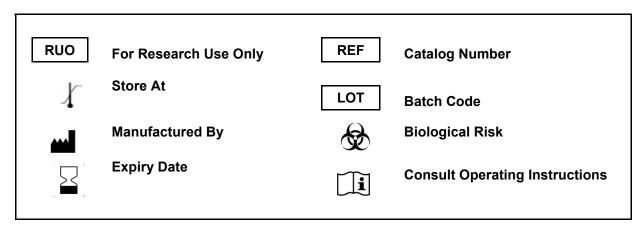

KRISHGEN BioSystems For US/Europe Customers: toll free +1(888)-970-0827 | tel +1(562)-568-5005 For Asia/India Customers: +91(22)-49198700 Email: sales@krishgen.com | http://www.krishgen.com

Human NT-Terminal Fragment of The Medullary Natriuretic Peptide, NT-ProBNP GENLISA™ ELISA



Introduction:

The GENLISA™ ELISA kits are used for assessing the specific biomarker in samples analytes which may be serum, plasma and cell culture supernatant as validated with the kit. The kit employs a sandwich ELISA technique which leads to a higher specificity and increased sensitivity compared to conventional competitive ELISA kits which employ only one antibody. Double antibodies are used in this kit.

Intended Use:

The Human NT-Terminal Fragment of The Medullary Natriuretic Peptide, NT-ProBNP GENLISA™ ELISA kit is used as an analytical tool for quantitative determination of Human NT-Terminal Fragment of The Medullary Natriuretic Peptide, NT-ProBNP in serum, plasma and other biological samples.

Principle:

The method employs sandwich ELISA technique. Monoclonal antibodies are pre-coated onto microwells. Samples and standards are pipetted into microwells and Human NT-Terminal Fragment of The Medullary Natriuretic Peptide, NT-ProBNP present in the sample are bound by the antibodies. Biotin labeled antibody is added and followed by Streptavidin:HRP is pipetted and incubated to form a complex. After washing microwells in order to remove any non-specific binding, the substrate solution (TMB) is added to microwells and color develops proportionally to the amount of Human NT-Terminal Fragment of The Medullary Natriuretic Peptide, NT-ProBNP in the sample. Color development is then stopped by addition of stop solution. Absorbance is measured at 450 nm.

Materials Provided:

- 1. Human NT-ProBNP Antibody Coated Microtiter Plate (8 x 12 wells) 1 no
- 2. Human NT-ProBNP Standard (lyophilized, concentrated, 5000 pg/ml) 2 vials
- 3. Biotinylated NT-ProBNP Antibody (concentrated) 120 ul
- 4. Streptavidin: HRP Conjugate (concentrated) 120 ul
- 5. Standard Diluent 20 ml
- 6. Biotin Antibody Dilution Buffer 12 ml
- 7. HRP Conjugate Dilution Buffer 12 ml
- 8. (20X) Wash Buffer 25 ml
- 9. TMB Substrate 12 ml
- 10. Stop Solution 12 ml
- 11. Instruction Manual

Materials to be provided by the End-User:

- 1. Microtiter Plate Reader able to measure absorbance at 450 nm.
- 2. Adjustable pipettes and multichannel pipettor to measure volumes ranging from 25 ul to 1000 ul
- 3. Deionized (DI) water
- 4. Wash bottle or automated microplate washer
- 5. Clean tubes and Eppendorf tubes
- 6. Precision single and multi-channel pipette and disposable tips.
- 7. 37°C incubator
- 8. Timer.

Handling/Storage:

- 1. All reagents should be stored as indicated on the component label.
- 2. All the reagents and wash solutions should be used within 12 months from manufacturing date.
- 3. Before using, bring all components to room temperature (18-25°C). Upon assay completion ensure all components of the kit are returned to appropriate storage conditions.
- 4. The Substrate is light-sensitive and should be protected from direct sunlight or UV sources.



Health Hazard Warnings:

- 1. Reagents that contain preservatives may be harmful if ingested, inhaled or absorbed through the skin.
- 2. For Research Use Only.



Sample Preparation and Storage:

Specimens should be clear and non-hemolyzed. Samples should be run at a number of dilutions to ensure accurate quantitation.

- 1. Extract as soon as possible after specimen collection as per relevant procedure. The samples should be tested as soon as possible after the extraction. Alternately the extracted samples can be kept in -20°C. Avoid repeated freeze-thaw cycles.
- 2. **Serum-** Coagulate at room temperature for 10-20 minutes; centrifuge for 20-min at 2000-3000 rpm. Remove the supernatant. If precipitation appears, recentrifuge.
- 3. **Plasma-** Use EDTA or citrate plasma as an anticoagulant, mix for 10-20 minutes; centrifuge for 15-min at 2000-3000 rpm. Remove the supernatant carefully. If precipitation appears, recentrifuge.
- 4. **Cell Culture Supernatant-** Collect sample in a sterile container. Centrifuge for 20-mins at 2000-3000 rpm. Remove the supernatant carefully. When examining the components within the cell, dilute cell suspension with PBS (pH 7.2-7.4), if cell concentration is greater than 1 million/ml. Damage the cells by repeated freeze-thaw cycles to release intracellular components. Centrifuge for 20-min at 2000-3000 rpm. If precipitation appears, centrifuge again.
- 5. **Tissue Samples-** Rinse tissues in PBS (pH 7.4) to remove excess blood thoroughly and weigh before homogenization. Mince tissues and homogenize them in PBS (pH7.4) with a glass homogenizer on ice. Thaw at 2-8°C or freeze at -20°C. Centrifuge at 2000-3000 RPM for approximately 20 minutes and collect the supernatant carefully.

Note: Grossly hemolyzed samples are not suitable for use in this assay.

Reagent Preparation (all reagents should be diluted immediately prior to use):

- 1. Label any aliquots made with the kit Lot No and Expiration date and store it at appropriate conditions mentioned.
- 2. Bring all reagents to Room temperature before use.
- 3. To make Wash Buffer (1X) 500 ml; dilute 25 ml of (20X) Wash Buffer in 475 ml of Dl water.
- 4. Streptavidin:HRP Conjugate & Biotinylated NT-ProBNP Antibody Working Solution Briefly spin or centrifuge the Streptavidin:HRP Conjugate & Biotinylated NT-ProBNP Antibody before use. Dilute them to the working concentration 100-fold with HRP Conjugate Dilution Buffer & Biotin Antibody Dilution Buffer, respectively.
- 5. **Standards Preparation**: Reconstitute original NT-Terminal Fragment of The Medullary Natriuretic Peptide, NT-ProBNP with 1.0 ml of Standard Diluent. Keep the standard for 10 mins with gentle agitation before making further dilutions. Prepare the additional Standards by serially diluting the standard stock solution as per the below table.

| Standard Concentration | Standard Vial | Dilution Particulars |
|------------------------|-------------------|--|
| 5000 pg/ml | Original Standard | Reconstitute with 1.0 ml Standard Diluent |
| 2500 pg/ml | Standard No.8 | 500 ul Original Standard + 500 ul Standard Diluent |
| 1250 pg/ml | Standard No.7 | 500 ul Standard No.8 + 500 ul Standard Diluent |
| 625 pg/ml | Standard No.6 | 500 ul Standard No.7 + 500 ul Standard Diluent |
| 312 pg/ml | Standard No.5 | 500 ul Standard No.6 + 500 ul Standard Diluent |
| 156 pg/ml | Standard No.4 | 500 ul Standard No.5 + 500 ul Standard Diluent |
| 78 pg/ml | Standard No.3 | 500 ul Standard No.4 + 500 ul Standard Diluent |
| 39 pg/ml | Standard No.2 | 500 ul Standard No.3 + 500 ul Standard Diluent |
| 0 pg/ml | Standard No.1 | 500 ul Standard Diluent only |

Procedural Notes:

Human NT-Terminal Fragment of The Medullary Natriuretic Peptide, NT-ProBNP GENLISA™ ELISA



- 1. In order to achieve good assay reproducibility and sensitivity, proper washing of the plates to remove excess un-reacted reagents is essential.
- 2. High Dose Hook Effect may be observed in samples with very high concentrations of Human NT-Terminal Fragment of The Medullary Natriuretic Peptide, NT-ProBNP. High Dose Hook Effect is due to excess of antibody for very high concentrations of Human NT-Terminal Fragment of The Medullary Natriuretic Peptide, NT-ProBNP present in the sample.
- 3. Human NT-Terminal Fragment of The Medullary Natriuretic Peptide, NT-ProBNP concentration of the undiluted sample is less than the diluted sample, this may be indicative of the Hook Effect.
- Avoid assay of Samples containing sodium azide (NaN₃), as it could destroy the HRP activity resulting in under-estimation of the amount of Human NT-Terminal Fragment of The Medullary Natriuretic Peptide, NT-ProBNP.
- 5. It is recommended that all Standards and Samples be assayed in duplicates or triplicates.
- 6. Maintain a repetitive timing sequence from well to well for all the steps to ensure that the incubation timings are same for each well.
- 7. If the Substrate has a distinct blue color prior to use it may have been contaminated and use of such substrate can lead to compromisation of the sensitivity of the assay.
- 8. The plates should be read within 30 minutes after adding the Stop Solution.
- 9. Make a work list in order to identify the location of Standards and Samples.

Assay Procedure:

- 1. It is strongly recommended that all Standards and Samples be run in duplicates or triplicates. A standard curve is required for each assay.
- 2. Add 100 ul prepared Standards and Samples to respective wells.
- 3. Cover the plate with a sealer and incubate for 60 minutes at 37°C.
- Aspirate the plate. DO NOT WASH.
- 5. Pipette 100 ul Biotinylated NT-ProBNP Antibody Working Solution to all wells.
- 6. Cover the plate with a sealer and incubate for 60 minutes at 37°C.
- 7. Aspirate and wash plate 4 times with diluted Wash Buffer (1X) and blot residual buffer by firmly tapping plate upside down on absorbent paper. Wipe of any liquid from the bottom outside of the microtiter wells as any residue can interfere in the reading step.
- 8. Pipette 100 ul Streptavidin:HRP Conjugate Working Solution to all wells. Mix well.
- 9. Cover the plate with a sealer and incubate for 30 minutes at 37°C.
- 10. Aspirate and wash as per Step (7) above.
- 11. Pipette 100 ul TMB Substrate in all the wells.
- 12. Incubate the plate at **37°C** for **10 minutes**. DO NOT SHAKE or else it may result in higher backgrounds and worse precision. Positive wells should turn bluish in color.
- 13. Pipette 100 ul of Stop Solution to all wells. The wells should turn from blue to yellow in color.
- 14. Read the absorbance at 450 nm with a microplate within 10-15 minutes after addition of Stop solution.

Calculation of Results:

Determine the Mean Absorbance for each set of duplicate or triplicate Standards and Samples. Using Graph Paper, plot the average value (absorbance 450nm) of each standard on the Y-axis versus the corresponding concentration of the standards on the X-axis. Draw the best fit curve through the standard points. To determine the unknown Human NT-Terminal Fragment of The Medullary Natriuretic Peptide, NT-ProBNP concentrations, find the unknown's Mean Absorbance value on the Y-axis and draw a horizontal line to the standard curve. At the point of intersection, draw a vertical line to the X-axis and read the Human NT-Terminal Fragment of The

Cat No#KBH4693, Ver3.0

Human NT-Terminal Fragment of The Medullary Natriuretic Peptide, NT-ProBNP GENLISA $^{\text{TM}}$ ELISA



Medullary Natriuretic Peptide, NT-ProBNP Concentration.

If samples were diluted, multiply by the appropriate dilution factor. Software which is able to generate a cubic spline curve-fit or 4-PL is best recommended for automated results

Note:

It is recommended to repeat the assay at a different dilution factor in the following cases:

- If the sample absorbance value is below the first standard.

Quality Control:

It is recommended that for each laboratory assay appropriate quality control samples in each run to be used to ensure that all reagents and procedures are correct.

Performance Characteristics of the Kit:

This kit has been validated. Please view the details herein below.

Standard Calibration Range:

39 pg/ml – 2500 pg/ml

Sensitivity:

Limit Of Quantification:

It is defined as the lowest detectable concentration that can be determined with an acceptable repeatability and the LOQ was found to be 17 pg/ml.

Specificity:

This assay has high sensitivity and excellent specificity for detection of NT-Terminal Fragment of The Medullary Natriuretic Peptide, NT-ProBNP. No significant cross-reactivity or interference between NT-Terminal Fragment of The Medullary Natriuretic Peptide, NT-ProBNP and analogues was observed.

Recovery

Matrices listed below were spiked with certain level of NT-Terminal Fragment of The Medullary Natriuretic Peptide, NT-ProBNP and the recovery rates were calculated by comparing the measured value to the expected amount of NT-Terminal Fragment of The Medullary Natriuretic Peptide, NT-ProBNP in samples.

| Matrix | Recovery Range (%) | Average (%) |
|---------------------|--------------------|-------------|
| serum(n=5) | 79-91 | 86 |
| EDTA plasma(n=5) | 90-103 | 97 |
| heparin plasma(n=5) | 85-95 | 90 |

Precision:

Intra-Assay: CV<10% Inter-Assay: CV<12%

Linearity

The linearity of the kit was assayed by testing samples spiked with appropriate concentration of NT-Terminal Fragment of The Medullary Natriuretic Peptide, NT-ProBNP and their serial dilutions. The results were demonstrated by percentage of calculated concentration to the expectation.

| Sample | 1:2 | 1:4 | 1:8 | 1:16 |
|---------------------|---------|---------|--------|---------|
| serum(n=5) | 88-97% | 82-103% | 80-96% | 79-101% |
| EDTA plasma(n=5) | 87-102% | 92-104% | 81-94% | 82-95% |
| heparin plasma(n=5) | 93-105% | 84-96% | 90-99% | 83-91% |

Safety Precautions:

- This kit is For Research Use only. Follow the working instructions carefully.
- The expiration dates stated on the kit are to be observed. The same relates to the stability stated for reagents
- Do not use or mix reagents from different lots.
- Do not use reagents from other manufacturers.

Human NT-Terminal Fragment of The Medullary Natriuretic Peptide, NT-ProBNP GENLISA™ ELISA



- · Avoid time shift during pipetting of reagents.
- All reagents should be kept in the original shipping container.
- Some of the reagents contain small amount of sodium azide (< 0.1 % w/w) as preservative. They must not be swallowed or allowed to come into contact with skin or mucosa.



- Source materials maybe derived from Human body fluids or organs used in the preparation of this kit were tested and found negative for HBsAg and HIV as well as for HCV antibodies. However, no known test guarantees the absence of such viral agents. Therefore, handle all components and all patient samples as if potentially hazardous.
- Since the kit contains potentially hazardous materials, the following precautions should be observed
- Do not smoke, eat or drink while handling kit material
- Always use protective gloves
- Never pipette material by mouth
- Wipe up spills promptly, washing the affected surface thoroughly with a decontaminant.
- In any case GLP should be applied with all general and individual regulations to the use of this kit.

Cat No#KBH4693, Ver3.0 www.krishgen.com 6



Typical Example of a Work List

| Well # | Contents | Absorbance at 450nm | Mean Absorbance | Interpolated Concentration |
|----------|---------------|---------------------|--------------------|----------------------------|
| 1A | Standard No.1 | | | |
| 2A | Standard No.1 | | | |
| 1B | Standard No.2 | | | |
| 2B | Standard No.2 | | | |
| 1C | Standard No.3 | | | |
| 2C | Standard No.3 | | | |
| 1D | Standard No.4 | | | |
| 2D | Standard No.4 | | | |
| 1E | Standard No.5 | | | |
| 2E | Standard No.5 | | | |
| 1F | Standard No.6 | | | |
| 2F | Standard No.6 | | | |
| 1G | Standard No.7 | | | |
| 2G | Standard No.7 | | | |
| 1H | Standard No.8 | | | |
| 2H | Standard No.8 | | | |
| 3A | Cample | | | |
| 4A | Sample | | | |
| 3B 4B | Sample | | | |

LIMITED WARRANTY

Krishgen Biosystems does not warrant against damages or defects arising in shipping or handling, or out of accident or improper or abnormal use of the Products; against defects in products or components not manufactured by Krishgen Biosystems, or against damages resulting from such non-Krishgen Biosystems made products or components. Krishgen Biosystems passes on to customer the warranty it received (if any) from the maker thereof of such non Krishgen made products or components. This warranty also does not apply to Products to which changes or modifications have been made or attempted by persons other than pursuant to written authorization by Krishgen Biosystems.

THIS WARRANTY IS EXCLUSIVE. The sole and exclusive obligation of Krishgen Biosystems shall be to repair or replace the defective Products in the manner and for the period provided above. Krishgen Biosystems shall not have any other obligation with respect to the Products or any part thereof, whether based on contract, tort, and strict liability or otherwise. Under no circumstances, whether based on this Limited Warranty or otherwise, shall Krishgen Biosystems be liable for incidental, special, or consequential damages.

This Limited Warranty states the entire obligation of Krishgen Biosystems with respect to the Products. If any part of this Limited Warranty is determined to be void or illegal, the remainder shall remain in full force and effect.

Krishgen Biosystems. 2023

THANK YOU FOR USING A KRISHGEN PRODUCT!

KRISHGEN BIOSYSTEMS®, GENLISA®, DHARMAPLEX™, GENBULK™, GENLISA™, KRISHZYME®, KRISHGEN®, KRIBIOLISA®, KRISHPLEX®, TITANIUM®, QUALICHEK® are registered trademarks of KRISHGEN BIOSYSTEMS. ©KRISHGEN BIOSYSTEMS. ALL RIGHTS RESERVED. KRISHGEN BIOSYSTEMS | OUR REAGENTS | YOUR RESEARCH |

Cat No#KBH4693, Ver3.0 www.krishgen.com 7



SYMBOLS KEY

| MTP | Coated Microtiter Plate (8x12 wells) |
|--------------|--------------------------------------|
| STD | Standard |
| BIOTIN AB | Biotinylated Antibody |
| HRP CONJ | Conjugate Horseradish Peroxidase |
| BIOTIN DIL | Biotin Antibody Dilution Buffer |
| HRP DIL | HRP Conjugate Dilution Buffer |
| STD DIL | Standard Diluent |
| 20X WASH BUF | (20X) Wash Buffer |
| SUB TMB | TMB Substrate |
| SOLN STOP | Stop Solution |
| []i | Consult Instructions for Use |
| REF | Catalog Number |
| | Expiration Date |
| * | Storage Temperature |

Cat No#KBH4693, Ver3.0 www.krishgen.com 8