KRIBIOLISA™ Anti-Eculizumab (SOLIRIS™) ELISA

REF : KBI2024

Ver 4.0

Enzyme Immunoassay for the Quantitative Determination of Anti-Eculizumab in human serum and plasma

| RUO | For Research Use Only | REF | Catalog Number |
|-----|-----------------------|----------|--------------------------------|
| X | Store At | LOT | Batch Code |
| 44 | Manufactured By | A | Biological Risk |
| | Expiry Date | Ĩ | Consult Operating Instructions |
| | | | |

For Research Purposes Only. Purchase does not include or carry the right to resell or transfer this product either as a stand-alone product or as a component of another product. Any use of this product other than the permitted use without the express written authorization of KRISHGEN BioSystems is strictly prohibited.





KRISHGEN BioSystems | For US / Europe: toll free +1(888)-970-0827 tel: +1(562)-568-500 For Asia / India: tel: +91(22)-49198700 Email: sales@krishgen.com

Introduction: Cat No#KBI2024 Ver 4.0


Eculizumab, sold under the trade name Soliris, is a medication used to treat paroxysmal nocturnal hemoglobinuria (PNH) and atypical hemolytic uremic syndrome (aHUS). It is a humanized monoclonal antibody functioning as a terminal complement inhibitor. In people with PNH, it reduces both the destruction of red blood cells and need for blood transfusion, but does not appear to affect the risk of death. ADA (Anrti-Drug Antibodies) has been observed in pre-clinical and clinical studies, resulting in significant changes in toxicology, pharmacokinetics and efficacy. These effects result from the generation of drug-induced (neutralizing) autoantibodies against Eculizumab and can be responsible for allergic reaction, or even anaphylactic shock.

Intended Use:

The KRIBIOLISA[™] Anti-Eculizumab (SOLIRIS[™]) ELISA is used as an analytical tool for quantitative determination of antibodies to Eculizumab in human serum and plasma.

Principle:

The method employs the quantitative sandwich enzyme immunoassay technique. Eculizumab is precoated onto microwells. Samples and standards are pipetted into microwells and antibodies to Eculizumab present in the sample are bound by the capture antibody. Then, a HRP (horseradish peroxidase) conjugated Eculizumab is pipetted and incubated. After washing microwells in order to remove any non-specific binding, the ready to use substrate solution (TMB) is added to microwells and color develops proportionally to the amount of Anti-Eculizumab in the sample. Color development is then stopped by addition of stop solution. Absorbance is measured at 450 nm.

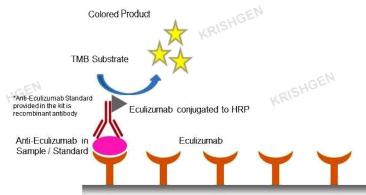

Materials Provided:

ELISA Coated Microplate

| Part | Description | Qty |
|------------------------------|--|-----------|
| Eculizumab Coated Microtiter | 96 well polystyrene microplate (12 strips of 8 wells) coated with | 1 x 96 |
| Plate | Eculizumab. | wells |
| | Recombinant Antibody to Eculizumab in a buffered protein base | |
| Anti-Eculizumab Standard | with protein stabilizer and preservatives 0.02% ethylisothiazolone and 0.02% bromonitrodioxane – lyophilized (1 ug/ml) | 2 vials |
| | Eculizumab conjugated to horseradish peroxidase with protein | |
| Eculizumab:HRP Conjugate | stabilizer and preservatives 0.02% methylisothiazolone and | 12 ml |
| | 0.02% bromonitrodioxane. | |
| Sample Diluent | Buffered protein base with protein stabilizer and preservatives | 2 x 50 ml |
| | 0.02% methylisothiazolone and 0.02% bromonitrodioxane | |
| | Buffered protein base with 1:1000 dilution human serum and | |
| Standard Diluent | protein stabilizer and preservatives 0.02% methylisothiazolone | 10 ml |
| | and 0.02% bromonitrodioxane | |
| | 25 ml/vial of a 20-fold concentrated solution of buffered | |
| (20X) Wash Buffer | surfactant with preservative thiomersol < 0.01%. May turn yellow | 25 ml |
| | over time. | |
| TMB Substrate | Stabilized Chromogen | 12 ml |
| Stop Solution | 2N Sulfuric Acid | 12 ml |
| Instruction Manual | | 1 no |

Materials to be provided by the End-User:

- 1. Microtiter Plate Reader able to measure absorbance at 450 nm.
- 2. Adjustable pipettes and multichannel pipettor to measure volumes ranging from 25 ul to 1000 ul
- 3. Deionized (DI) water
- 4. Wash bottle or automated microplate washer
- 5. Standard (mm/mm) graph paper or software for data analysis
- 6. Timer
- 7. Absorbent Paper

Handling/Storage:

- 1. All reagents should be stored at 2°C to 8°C for stability.
- 2. All the reagents and wash solutions should be used within 12 months from manufacturing date.
- 3. Before using, bring all components to room temperature (18-25°C). Upon assay completion ensure all components of the kit are returned to appropriate storage conditions.
- 4. The Substrate is light-sensitive and should be protected from direct sunlight or UV sources.

Health Hazard Warnings:

- 1. Reagents that contain preservatives may be harmful if ingested, inhaled or absorbed through the skin.
- 2. For Research use only.

Sample Preparation and Storage:

Blood is taken by venipuncture. Serum is separated after clotting by centrifugation. Plasma can be used, too. Lipaemic, hemolytic or contaminated samples should not be run. Repeated freezing and thawing should be avoided. If samples are to be used for several assays, initially aliquot samples and keep at - 20°C.

For Cell Culture Supernatant - If necessary, centrifuge to remove debris prior to analysis. Samples can be stored at -20°C or -80°C. Avoid repeated freeze-thaw cycles.

Preparation Before Use:

Allow samples to reach room temperature prior to assay. Take care to agitate patient samples gently in order to ensure homogeneity.

Serum and Plasma Test Sample preparation - Samples have to be diluted 1:1000 (v/v), e.g. for 1:1000 (1 ul sample + 999 ul **Sample Diluent**) prior to assay. The samples may be kept at 2 - 8°C for up to three days. Long-term storage requires -20°C.

Reagent Preparation (all reagents should be diluted immediately prior to use):

- 1. Label any aliquots made with the kit Lot No and Expiration date and store it at appropriate conditions mentioned.
- 2. Bring all reagents to Room temperature before use.
- 3. To make Wash Buffer (1X); dilute 25 ml of 20X Wash Buffer in 475 ml of DI water.
- 4. Standards Preparation: Reconstitute the concentrated Standard lyophilized vial with 1 ml of Standard Diluent to obtain a concentration of 1 ug/ml. Keep the vial for 15 mins with gentle agitation before making further dilutions. Dilute 320 ul of original Standard (1 ug/ml) with 180 ul of Standard Diluent to generate a 640 ng/ml Standard Solution. Prepare further Standards by serially diluting the Standard Solution as per the below table. Use the Standard Diluent as the Zero Standard (Standard No.0).

| Concentration | | |
|---------------|-------------------------|---|
| 1 ug/ml | Lyophilized Standard | Lyophilized Standard provided in the Kit + 1 ml of Standard Diluent |
| 640 ng/ml | Standard No.6 | 320 ul Reconstituted Standard (1 ug/ml) + 180 ul Standard Diluent |
| 320 ng/ml | Standard No.5 | 250 ul Standard No.6 + 250 ul Standard Diluent |
| 160 ng/ml | Standard No.4 | 250 ul Standard No.5 + 250 ul Standard Diluent |
| 80 ng/ml | Standard No.3 | 250 ul Standard No.4 + 250 ul Standard Diluent |
| 40 ng/ml | Standard No.2 | 250 ul Standard No.3 + 250 ul Standard Diluent |
| 20 ng/ml | Standard No.1 | 250 ul Standard No. 2 + 250 ul Standard Diluent |
| 0 ng/ml | Standard No.0 | Only Standard Diluent |
| <u> </u> | | |

Use the Standards as soon as possible upon reconstitution. Discard balance standard after use.

Procedural Notes:

- 1. In order to achieve good assay reproducibility and sensitivity, proper washing of the plates to remove excess un-reacted reagents is essential.
- 2. High Dose Hook Effect may be observed in samples with very high concentrations of Anti-Eculizumab. High Dose Hook Effect is due to excess of antibody for very high concentrations of Anti-Eculizumab present in the sample. High Dose Hook effect is most likely encountered from samples early in the purification process. If Hook Effect is possible, the samples to be assayed should be diluted with a compatible diluent. Thus if the Anti- Eculizumab concentration of the undiluted sample is less than the diluted sample, this may be indicative of the Hook Effect.
- 3. Avoid assay of Samples containing sodium azide (NaN₃), as it could destroy the HRP activity resulting in under-estimation of the amount of Anti- Eculizumab.
- 4. It is recommended that all Standards and Samples be assayed in duplicates.
- 5. Maintain a repetitive timing sequence from well to well for all the steps to ensure that the incubation timings are same for each well.
- 6. If the Substrate has a distinct blue color prior to use it may have been contaminated and use of such substrate can lead to compromisation of the sensitivity of the assay.
- 7. The plates should be read within 30 minutes after adding the Stop Solution.
- 8. Make a work list in order to identify the location of Standards and Samples.

Assay Procedure:

- 1. It is strongly recommended that all Standards and Samples be run in duplicates. A standard curve is required for each assay. All steps must be performed at 37°C
- 2. Pipette 100 ul of prepared Standards or diluted Samples into the respective wells
- 3. Cover the plate and incubate for 60 minutes at 37°C
- 4. Aspirate and wash plate 4 times with **Wash Buffer (1X)** and blot residual buffer by firmly tapping plate upside down on absorbent paper. Wipe of any liquid from the bottom outside of the microtiter wells as any residue can interfere in the reading step.
- 5. Pipette 100 ul of Eculizumab: HRP Conjugate into the respective wells.
- 6. Cover the plate and incubate for 60 minutes at 37°C.
- 7. Aspirate and wash plate 4 times with **Wash Buffer (1X)** same as in step 4.
- 8. Add 100 ul of TMB Substrate in each well.
- 9. Incubate the plate at 37°C for 30 minutes in dark. DO NOT SHAKE or else it may result in higher backgrounds and worse precision. Positive wells should turn bluish in color.
- 10. Pipette out **100 ul** of **Stop Solution**. Wells should turn from blue to yellow in color.
- 11. Read the absorbance at 450 nm with a microplate reader.

Calculation of Results:

Determine the Mean Absorbance for each set of duplicate or triplicate Standards and Samples. Using Semi-Log graph paper, plot the average value (absorbance 450 nm) of each standard on the Y-axis versus the corresponding concentration of the standards on the X-axis. Draw the best fit curve through the standard points. To determine the unknown Anti-Eculizumab concentrations, find the unknown's Mean Absorbance value on the Y-axis and draw a horizontal line to the standard curve. At the point of intersection, draw a vertical line to the X-axis and read the Anti-Eculizumab Concentration. If samples were diluted, multiply by the appropriate dilution factor.

Software which is able to generate a cubic spline curve-fit or a polynomial curve (2nd order) is best recommended for automated results.


Note:

It is recommended to repeat the assay at a different dilution factor in the following cases:

- If the sample absorbance value is below the first standard.
- If the absorbance value is equivalent or higher than the 640 ng/ml standard.

| Standard Concentration (ng/ml) | Absorbance A | Absorbance B | Mean Absorbance | Interpolated Concentration | % Interpolated Concentration against Actual Concentration |
|--------------------------------------|-----------------|-----------------|--------------------|-------------------------------|--|
| 0 | 0.091 | 0.100 | 0.095 | | |
| 20 | 0.112 | 0.127 | 0.120 | 19.9 | 99.6 |
| 40 | 0.137 | 0.147 | 0.142 | 41.3 | 103.2 |
| 80 | 0.243 | 0.275 | 0.259 | 90.8 | 113.5 |
| 160 | 0.494 | 0.577 | 0.536 | 155.9 | 97.4 |
| 320 | 1.503 | 1.472 | 1.488 | 321.1 | 100.4 |
| 640 | 2.776 | 2.669 | 2.723 | 639.5 | 99.9 |

Typical Data

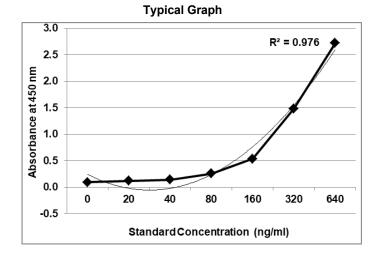

Quality Control:

It is recommended that for each laboratory assay appropriate quality control samples in each run to be used to ensure that all reagents and procedures are correct.

Performance Characteristics of the Kit:

This kit has been validated as per EMA/FDA guidelines in line with ICH Code for Harmonization of Biological Assays.

Sensitivity:

Limit Of Detection: It is defined as the lowest detectable concentration corresponding to a signal of Mean of '0' standard plus 2* SD.

10 replicates of '0' standards were evaluated and the LOD was found to be less than 18 ng/ml

Specificity:

The antibodies used in the kit are monoclonal antibodies, anti-idiotypic and specific for Anti-Eculizumab. The calibrators / standards used are calibrated against commercially source.

Precision:

Precision is defined as the percent coefficient of variation (%CV) i.e. standard deviation divided by the mean and multiplied by 100. Assay precision was determined by both intra (n=5 assays) and inter assay (n=5 assays) reproducibility on two pools with low (20 ng/ml), medium (80 ng/ml) and high (640 ng/ml) concentrations. While actual precision may vary from laboratory to laboratory and technician to technician, it is recommended that all operators achieve precision below these design goals before reporting results.

| Pool | Intra Assay %CV | Inter Assay %CV |
|--------|-----------------|-----------------|
| Low | <10% | <10% |
| Medium | <5% | <5% |
| High | <5% | <5% |

Safety Precautions:

- This kit is for Research Use only. Follow the working instructions carefully.
- The expiration dates stated on the kit are to be observed. The same relates to the stability stated for reagents
- Do not use or mix reagents from different lots.
- Do not use reagents from other manufacturers.
- · Avoid time shift during pipetting of reagents.
- All reagents should be kept in the original shipping container.
- Some of the reagents contain small amount of sodium azide (< 0.1 % w/w) as preservative. They must not be swallowed or allowed to come into contact with skin or mucosa.
- Source materials maybe derived from human body fluids or organs used in the preparation of this kit were tested and found negative for HBsAg and HIV as well as for HCV antibodies. However, no known test guarantees the absence of such viral agents. Therefore, handle all components and all patient samples as if potentially hazardous.
- Since the kit contains potentially hazardous materials, the following precautions should be observed
- Do not smoke, eat or drink while handling kit material
- Always use protective gloves
- Never pipette material by mouth
- Wipe up spills promptly, washing the affected surface thoroughly with a decontaminant.
- In any case GLP should be applied with all general and individual regulations to the use of this kit.

References:

Immunogenicity of meningococcus C vaccination in a patient with atypical hemolytic uremic syndrome (aHUS) on eculizumab therapy... Zlamy M, Hofer J, Elias J, Vogel U, Frosch M, Jungraithmayr T, Zimmerhackl LB, Prelog M...Pediatr Transplant... 2012...Wiley


Stratification of responders towards eculizumab using a structural epitope mapping strategy...Anna-Luisa Volk, Francis Jingxin Hu, Magnus M. Berglund, Erik Nordling, Patrik Strömberg, Mathias Uhlen & Johan Rockberg...Scientific reports...2016...Nature

Therapeutic drug monitoring of eculizumab: Rationale for an individualized dosing schedule...Philippe Gatault,Guillaume Brachet,David Ternant, Danielle Degenne, Guillaume Récipon, Christelle Barbet, Emmanuel Gyan, Valérie Gouilleux-Gruart, Cécile Bordes, Alexandra Farrell, Jean Michel Halimi, and Hervé Watier...MAbs...2015...Taylor & Francis

Sensitive, reliable and easy-performed laboratory monitoring of eculizumab therapy in atypical hemolytic uremic syndrome...Elena B.,VolokhinaaNicole C.A.,J.van de KaraGrethe,Bergsethb,Thea J., A.M.van der Veldena,Dineke Westraa ,Jack F.M.Wetzelsc,Lambertus P.van den Heuvelade, Tom EirikMollnes...Science direct...2015...Elsevier

Development and validation of an enzyme-linked immunosorbent assay to measure free eculizumab concentration in serum... Christophe Passot, Céline Desvignes, David Ternant, Theodora Bejan-Angoulvant, Anne-Claire Duveau, Philippe Gatault & Gilles Paintaud... Bioanalysis...2017...Future science

SCHEMATIC ASSAY PROCEDURE

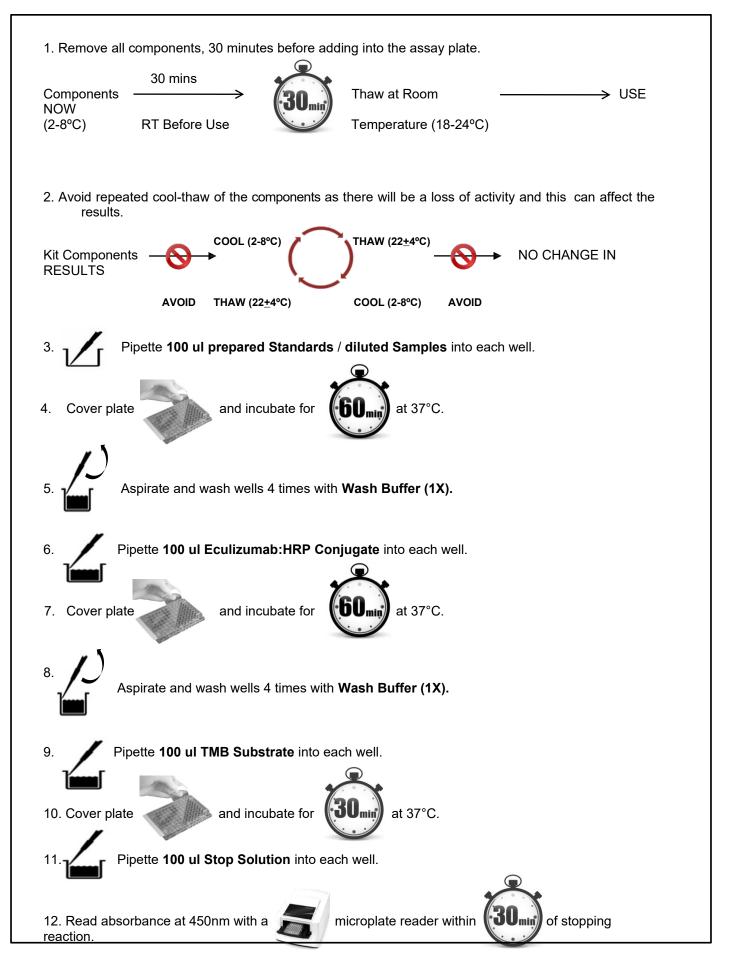

| Well # | Contents | Absorbance at 450nm | Mean Absorbance | ng/ml Anti- Eculizumab equivalent |
|--------|-----------|------------------------|--------------------|---|
| 1A | zero std | | | |
| 2A | zero std | | | |
| 1B | 20 ng/ml | | | |
| 2B | 20 ng/ml | | | |
| 1C | 40 ng/ml | | | |
| 2C | 40 ng/ml | | | |
| 1D | 80 ng/ml | | | |
| 2D | 80 ng/ml | | | |
| 1E | 160 ng/ml | | | |
| 2E | 160 ng/ml | | | |
| 1F | 320 ng/ml | | | |
| 2F | 320 ng/ml | | | |
| 1G | 640 ng/ml | | | |
| 2G | 640 ng/ml | | | |
| 1H | Sample | | | |
| 2H | Sample | | | |

Typical Example of a Work list

LIMITED WARRANTY

Krishgen Biosystems does not warrant against damages or defects arising in shipping or handling, or out of accident or improper or abnormal use of the Products; against defects in products or components not manufactured by Krishgen Biosystems, or against damages resulting from such non-Krishgen Biosystems made products or components. Krishgen Biosystems passes on to customer the warranty it received (if any) from the maker thereof of such non Krishgen made products or components. This warranty also does not apply to Products to which changes or modifications have been made or attempted by persons other than pursuant to written authorization by Krishgen Biosystems.

THIS WARRANTY IS EXCLUSIVE. The sole and exclusive obligation of Krishgen Biosystems shall be to repair or replace the defective Products in the manner and for the period provided above. Krishgen Biosystems shall not have any other obligation with respect to the Products or any part thereof, whether based on contract, tort, and strict liability or otherwise. Under no circumstances, whether based on this Limited Warranty or otherwise, shall Krishgen Biosystems be liable for incidental, special, or consequential damages.

This Limited Warranty states the entire obligation of Krishgen Biosystems with respect to the Products. If any part of this Limited Warranty is determined to be void or illegal, the remainder shall remain in full force and effect.

Krishgen Biosystems. 2023

THANK YOU FOR USING KRISHGEN PRODUCT!

KRISHGEN BIOSYSTEMS®, GENLISA®, DHARMAPLEX™, GENBULK™, GENLISA™, KRISHZYME®, KRISHGEN®, KRIBIOLISA®, KRISHPLEX®, TITANIUM®, QUALICHEK® are registered trademarks of KRISHGEN BIOSYSTEMS. ©KRISHGEN BIOSYSTEMS. ALL RIGHTS RESERVED.

KRISHGEN BIOSYSTEMS | OUR REAGENTS | YOUR RESEARCH |

Soliris[™] is the registered trade mark of Astra Zeneca

SYMBOLS KEY

| SHRPMICPNJ 283 NEWEDODF |
|----------------------------|
|----------------------------|

| | Eculizumab Coated Microtiter Plate (12x8 wells) |
|----|---|
| | Anti-Eculizumab Standard lyophilized |
| | Conjugate Horseradish Peroxidase |
| | Sample Diluent |
| | Standard Diluent |
| | (20X) Wash Buffer |
| | TMB Substrate |
| | Stop Solution |
| i | Consult Instructions for Use |
| | Catalog Number |
| | Expiration Date |
| ľ. | Storage Temperature |